

## **Practicum 4, Fall 2010**

### **The acquisition time: weighting and linear prediction**

Strychnine, dissolved in CDCl<sub>3</sub>

This demonstration develops the inverse relationship between time and frequency as we go back and forth between the FID, which occurs in time (it is time-domain data) and the NMR spectrum, which is displayed over a frequency axis (it is frequency-domain data).

Weighting functions are developed as a way to decrease the contribution to the spectrum of noise, and to correct for data truncation. They can also be used to enhance resolution, provided the data are of good quality.

### **Improving signal/noise with weighting**

Again, we begin by acquiring a nice 1-dimensional spectrum.

Now click on the FID display button [1] to look at the corresponding FID (**Figure 1**).

Note that the axis is now a time axis [2], and it runs from zero to the value you used as the acquisition time (at) [3].

For best spectral resolution, the FID should have decayed to zero by the time at ends. However making at much longer collects noise only, and decreases the quality of the spectrum, because every point in the FID is combined to generate each point in the spectrum, by the Fourier transformation. Look at the X-axis while displaying the FID. Decide at what time the signal is essentially gone, and choose that time for use as your improved acquisition time [4].

Thus, in this instance, a better choice of at is 1.6 seconds.

You can make this decision at the same time as you adjust the sweep width, based on a trial spectrum.

You can also cause the software to simply not use all of the digitized data as input to the FID, and you can do so in a tapered way.

Under Processing>Weighting (**Figure 2**, [0]) activate 'Interactive Weighting' [1]. In the graphical display window, the bottom panel shows the FID [2], the middle panel shows the function that will be multiplied point-for-point with the FID [3] and the top panel shows what the Fourier transform of this product will be [4]. Neither of the latter two panels display anything yet.

Click next to 'sinebell' to activate a sine function [5].

Now the middle panel shows a sine function (**Figure 3**, [1]). Click in the upper panel using the RightMouse (RM) [2]. The anticipated spectrum appears as a preview (not shown).

The spectrum that results has distorted peak shapes (**Figure 4 [1]**). Go to Process>Display, click 'Transform' and 'Display Spectrum' to get a better view. Click on 'Transform' in order to perform a new weighted Fourier transform incorporating the weighting function you just chose (**Figure 5**).

The strong negative excursions on each side of each signal are because the first point of the weighted FID will be zero, because the first point of the weighting function is zero. The value of the first point in the time domain determines the total integral under the spectrum after Fourier transformation. Thus in the sine-weighted spectrum, all the peaks have as much negative area as positive area, for a total area of zero, in theory.

To prevent the first point from having a zero value while nonetheless exploiting the sine function's decay to zero at long times, we shift it. In Processing>Weighting activate 'Interactive Weighting' [1] again and activate 'shift' [2] (**Figure 6**). Activate the spectrum preview by clicking in the top panel with the RM [3]. For shift, enter -1 \* the value you have next to 'sinebell'. Now you see that the peaks are broader but no longer have negative dips on both sides. Their integrals are positive. To shift the sine bell so that its peak is at time=0, set the shift equal to -1 \* the sine width.

The lines are broadened now because data are being squeezed to zero by the weighting function at time equal to the value used for the 'sinebell' parameter (look at the time axis in the lower panel of the interactive weighting display [4]).

After processing with your new weighting function, the 'shifted sinebell', determine the width at half-height of the signals in the frequency-domain spectrum by going to Process>Cursors/LineLists, and clicking on 'Transform' (**Figure 7**). Use the tools on the RHS to zoom in on a representative peak [0], place the cursor on it and click 'Place on nearest line' [1] then 'Show linewidth' [2]. In the space above the command line, as well as the RH response box at the bottom, the linewidth at half height is provided (2.55 Hz) along with the digital resolution [3].

Note that this screen also gives you the possibility of making the chosen line the center of the next spectrum you collect, by clicking on 'Move transmitter' [4]. This will change the value of tof (carrier frequency = 'transmitter offset') without changing the sweep width sw.

Go back to Processing>Weighting and activate 'interactive Weighting' again. Now set the shift to only -1\* *half* the value in sinebell. This is a 45 ° shifted sine bell, as opposed to the previous 90 ° shifted sinebell. The peaks will now be a little shorter, but their linewidths should be a little longer. This weighting function also has the beneficial effect of cancelling out 'feet' that often accompany peaks, complicating integrals and causing streaks in 2D spectra. Thus, it is often used.

You can also select other weighting functions. The Gaussian function and the line broadening (an exponential decrease) both have full amplitude at the first point and therefore do not require shifting, although shifting is often helpful (see above) and is possible with the Gaussian. Try the different options !

### Correcting for short $\Delta t$ with weighting

Sometimes we intentionally do not digitize the FID for as long as it takes to decay to zero. Sometimes this is because the sample contains one species with a very long  $T_2$  and that species is not of interest to us so we have optimized the experiment parameter in favour of other species present. Sometimes it is to keep the data set size smaller and easier to manipulate (eg. for 2D spectra). The result of having stopped data collection before the FID is really over is shown in **Figure 8**. In this case, the  $\Delta t$  chosen was longer than the  $T_2$  of the major constituent of the sample (strychnine), but it was shorter than the  $T_2$  of the TMS. Thus we see 'sinc wiggles' (feet) at the base of the TMS line [1] which is sharp (long  $T_2$ ). These are a 'truncation artifact'. (Note that the lines of strychnine which are broader (shorter  $T_2$ ) do not have sinc wiggles [2].)

There are two ways to 'fix' truncation artifacts without re-collecting the data. One is to apply a weighting function that compresses the data to zero at the time when the data collection ceased. The other is to apply linear prediction to extend the acquired data with synthetic 'predicted' data.

**Figure 9** shows the Process>Weighting interactive display with the truncated FID. (Note how the oscillations continue at the end of the data.)

Activate a Gaussian weighting function (**Figure 10**, [1]), click in the middle panel using the LM to extend or compress the Gaussian interactively [2] (the horizontal position of the mouse matters but not the vertical position within the middle panel) and note the effect on the preview of the frequency domain spectrum in the top panel [3] of the interactive weighting display. The sinc wiggles are gone, but the sharp line is broadened. In each instance, you will have to decide which is more important: no sinc wiggles or sharper lines.

Click on 'Transform' to apply your new weighting function and Fourier transform the spectrum. (**Figure 11** [1]). Compare the result with the spectrum in figure 8 that was not weighted prior to transformation.

### Linear Prediction

If the  $\Delta t$  is shorter than the decay of the data, you can have sinc wiggles. If you want to remove those without broadening your lines, you need more time-domain data (narrower freq-domain lines need longer time-domain data).

The software 'knows' that the data you collected should be the sum of many cosine functions each with the frequency of one peak and each damped according to that peak's  $T_2$  ( $1/\pi\Delta\nu$ ). Thus, if you tell it how many strong lines there are in your spectrum (approximately), the software will fit your acquired data in terms of a sum of damped Cos' to obtain sets of  $\nu$  and  $T_2$  s. (It will also generate phases and amplitudes for each.) Then the software uses this set of parameters to calculate what the 'ideal' FID should be for any time points you want. If we want to extend a short FID, we will ask for 'forward' linear prediction and we will ask the software to predict the perfect FID for time points after we stopped actually digitizing data.

In **Figure 12**, the actual data are in black and the predicted data are in blue.

To generate synthetic (predicted) points that double the data, go to Process>More 1D (**Figure 13**).

A number of decisions are required. We need to tell the software whether to linear predict forward in time or backward (parameter lpop = 'f' or 'b') (**Figure 13** [1]).

We need to specify how many peaks are to be assumed in fitting our data and which actual data points are to be the basis for the fit to determine peak  $\nu$  and  $T_2$  for each peak. Having seen your spectrum, estimate the number of *strong* peaks, for example 16 (don't be too greedy). Enter this next to 'coeff' [2] (parameter = lpfilt). For forward linear prediction we do this by specifying the number of points to be used [3] (parameter lpnupts can be as large as np, the number of points you have), and the *last* data point to be used (strtlp, 'start lp' in Figure 13 [4]). Thus if you want to use points 1-609, give 609 as the number of basis points [3] and 609 as the last point [4]. (This is called the 'starting point', but the software looks *backward* from that point when doing forward linear prediction. Just one of those things . . .)

If you chose to simulate 'n' signals you will need a minimum of  $2*n$  data points, but you should really use at least twice that many. I chose to generate points based on a model dominated by 16 signals. This requires that I use at least 32 points as a basis for my simulation. In our case, with nice strong data, we can afford to use many more points than that. To use the entire first half of the FID, use 609 points, from point 609 back to the start.

Enter the number of 'predicted points' you would like to add to your FID [5] (parameter lpext). To double the length of your FID enter the same number as you see below in the box 'Acquired points' [6]. This number is fixed at the time you collect data, it is reported here for your convenience (parameter name np). If these are to be added after the last actual data point, then the first predicted point is  $1 + np$ . Enter this value next to 'starting at' [7] (parameter strtext).

There is another choice related to your Fourier transformation operation. Even though we have  $np=1218$ , we can choose any power of two for the output of our Fourier transform. To use all 1218 input points, we would choose to output to 2048 (next biggest power of 2). What is done is that the input is augmented to 2048 points too by addition of a string of zeros after the last data point (zero-filling). When we linear predict, we will be doing a better job, but putting simulated data there instead. However if our simulated points make the FID so much longer that it now exceeds 2048 points, we may increase the Fourier transform size to 4096, or 4k [8] (parameter fn, or Fourier number).

This will determine the digital resolution, the number of Hz per point. Since you would like digital resolution better or equal to your spectral resolution (that is to say a smaller number of Hz / point), make fn large enough that the digital resolution is a smaller number than  $1/2$  your sharpest line's width at half height. These values can be checked on a line by going to Process>Weighting, placing the cursor on the

line and clicking on 'Display Linewidth' (**Figure 11** [4]). (Alternately, type **dres**). You don't need to collect new data if **fn** is too small, just increase it and Fourier transform again by clicking on 'Transform'.


**Figure 14** shows the spectrum that was acquired with a very short 0.3 s at, but augmented with linear prediction to double the FID length. In this case no weighting was applied in order to show the effect of linear prediction alone. The analogous spectrum that was not linear predicted is shown in **Figure 15**.

To get the best possible spectrum from the short-at data set, combine the use of linear prediction and weighting, by activating weighting in the linear prediction window (**Figure 15** [1]). (go to Processing>Weighting to double the Gaussian length in accordance with your new doubled FID length. If you used a shifted Gaussian, also double the shift.) You can double these in the command line by typing **gf=2\*gf** (Gaussian function is now twice the previous Gaussian function) and **gfs=2\*gfs** (Gaussian shift is now twice the previous Gaussian shift).

In figures 13, 14 and 15 compare both the line widths of the signals and the magnitude of truncation artifacts (sinc wiggles at the base of sharp signals).


For completeness' sake, **Figure 16** shows how the meanings of the linear prediction parameters change when one uses linear prediction to fix early points, by 'backward' linear prediction (black data are actual data to be used as the basis for a fit and blue are data that are predicted and will overwrite bad initial data points). In this case the choice of the number of signals to model ('coef' or **lpfilt**) and number of actual data are the same, but now the data points to be used as the basis for the fit extend from the first point to be used ('starting at' = **strtlp**) out to later times by a number of points called 'basis pts' or **lpnpts**. The first predicted point will be at 'starting at' (**strtext**) and will extend backwards for a number of points you enter as 'predicted points' (**lpext**)

Figure 1  
Weighting & Linear  
Prediction

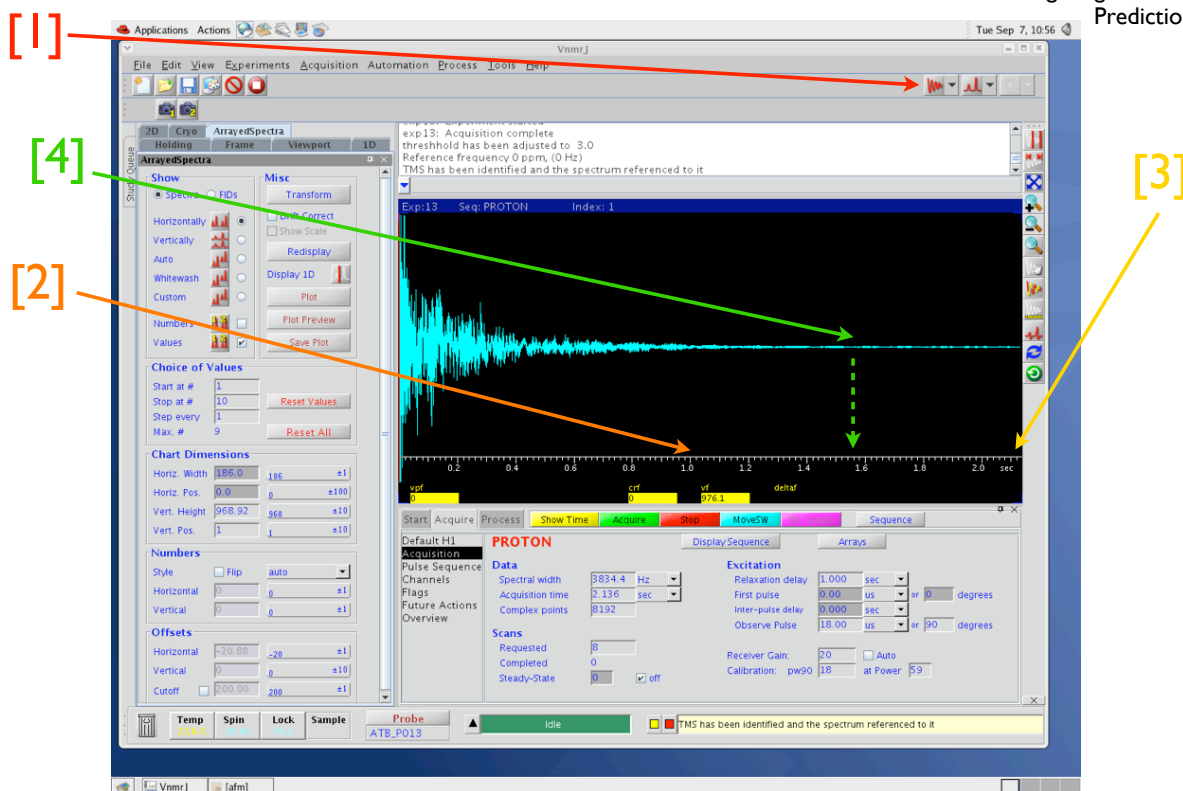


© A.-F. Miller 2010

Figure 2  
Weighting & Linear  
Prediction

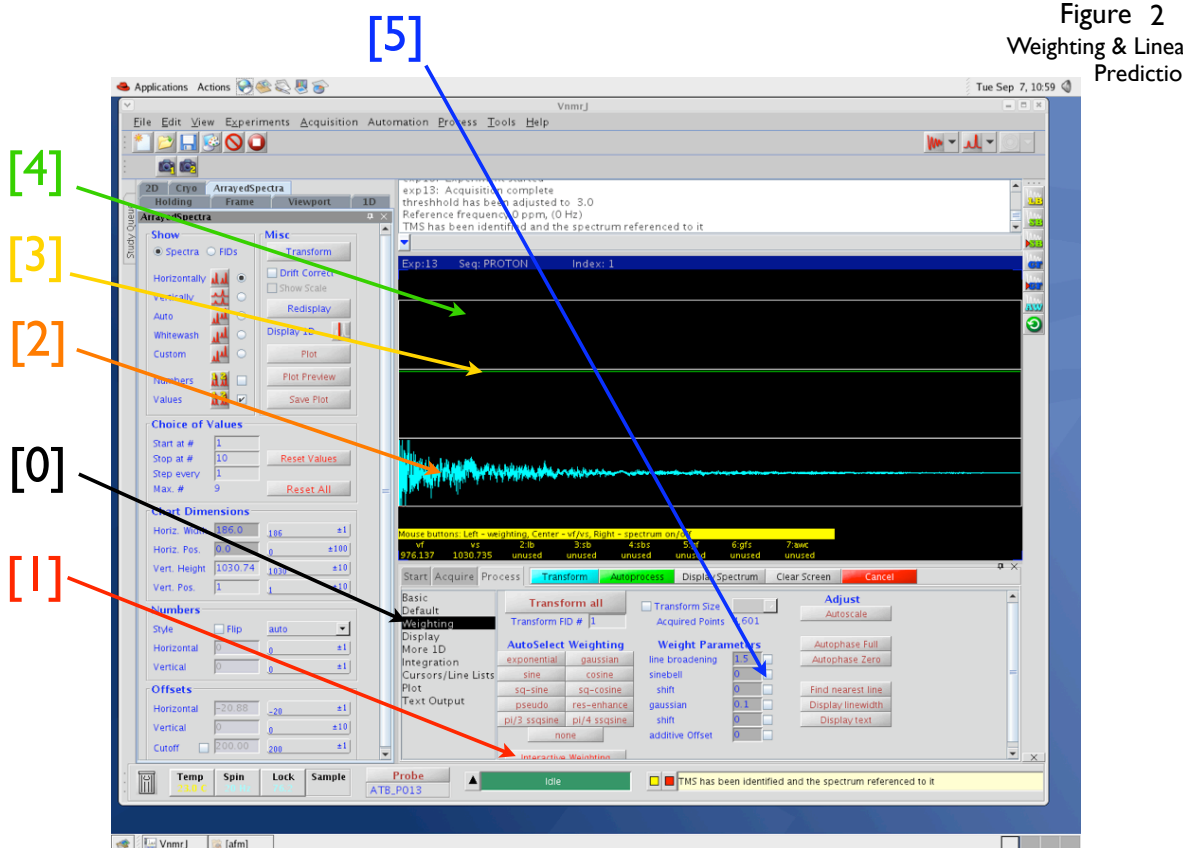


© A.-F. Miller 2010

Figure 3  
Weighting & Linear  
Prediction

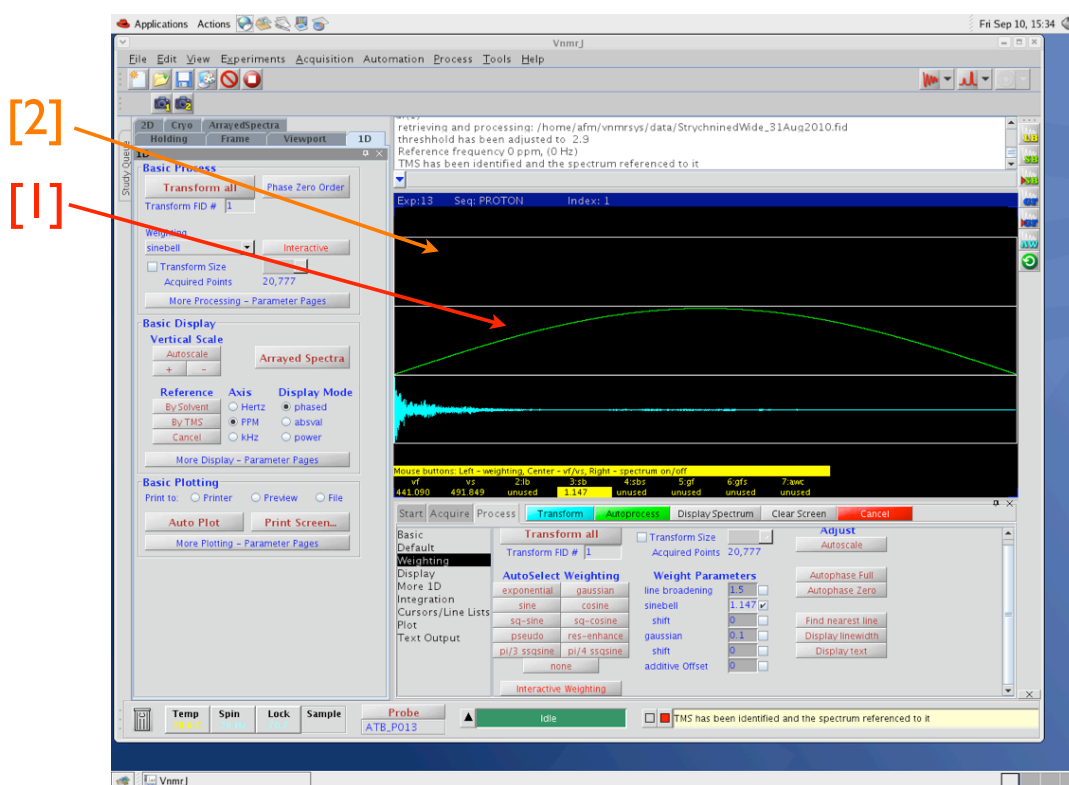


© A.-F. Miller 2010

Figure 4  
Weighting & Linear  
Prediction

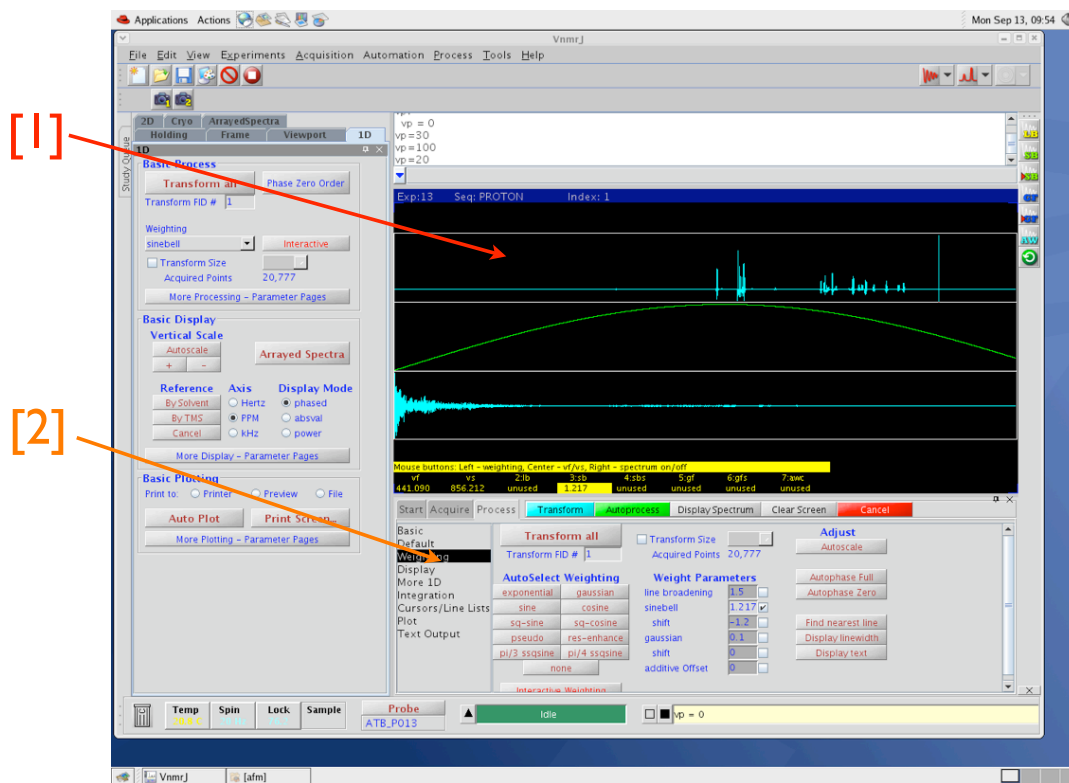


© A.-F. Miller 2010

Figure 5  
Weighting & Linear  
Prediction

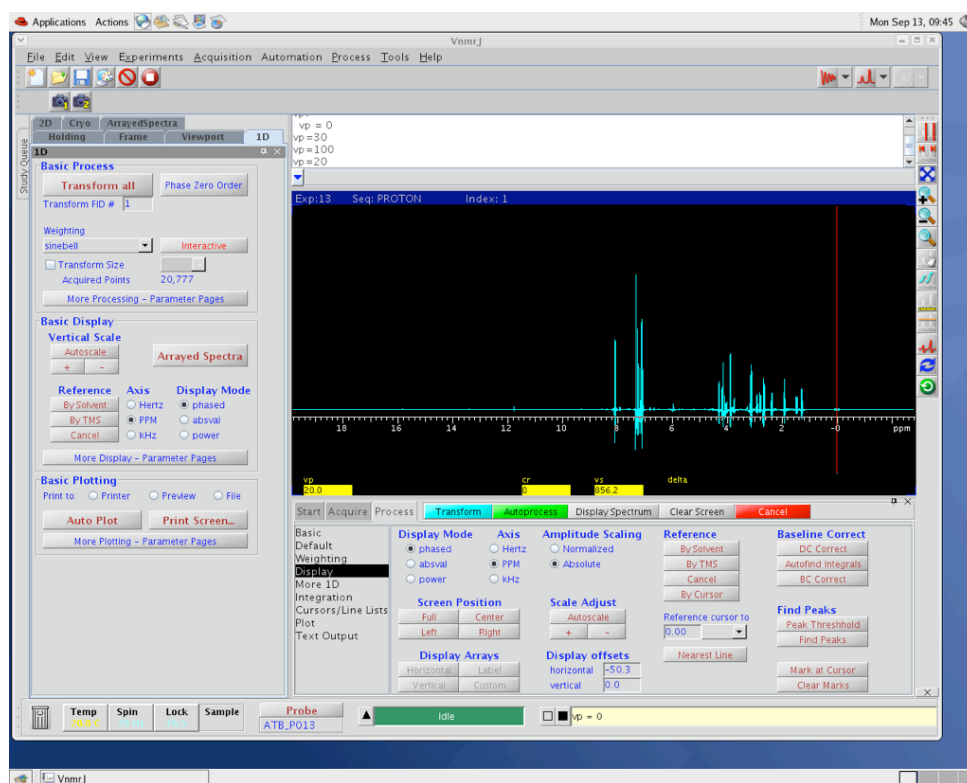

© A.-F. Miller 2010


Figure 6  
Weighting & Linear  
Prediction

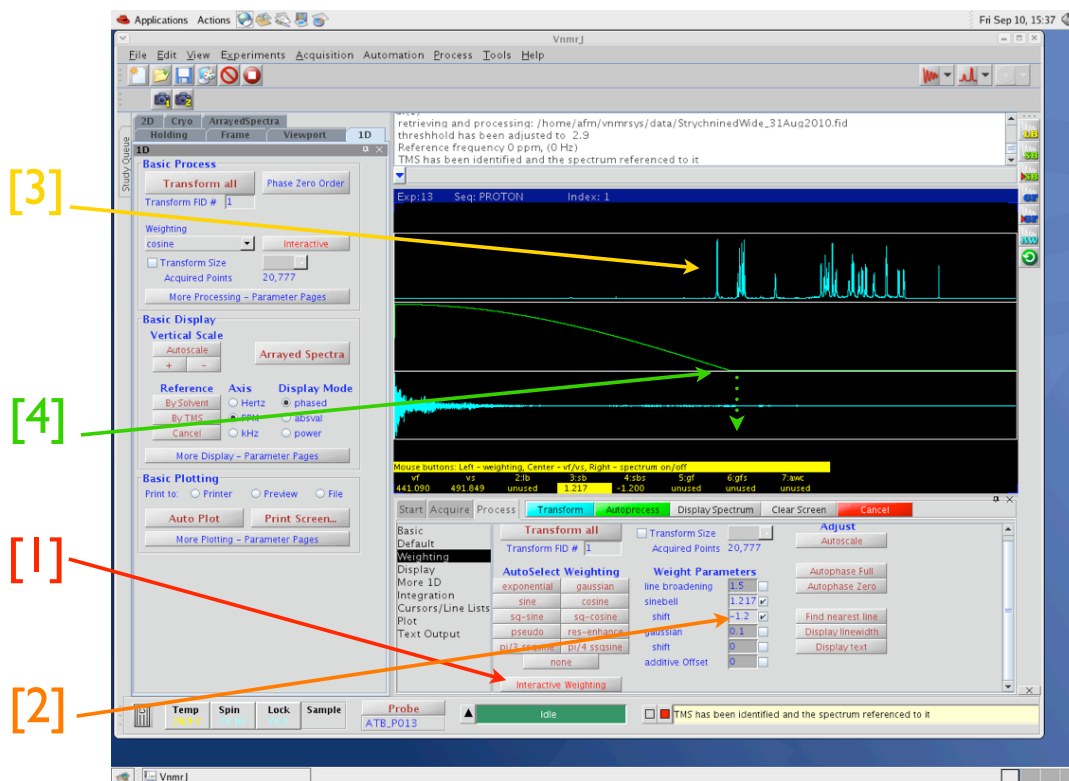

© A.-F. Miller 2010




Figure 7  
Weighting & Linear Prediction

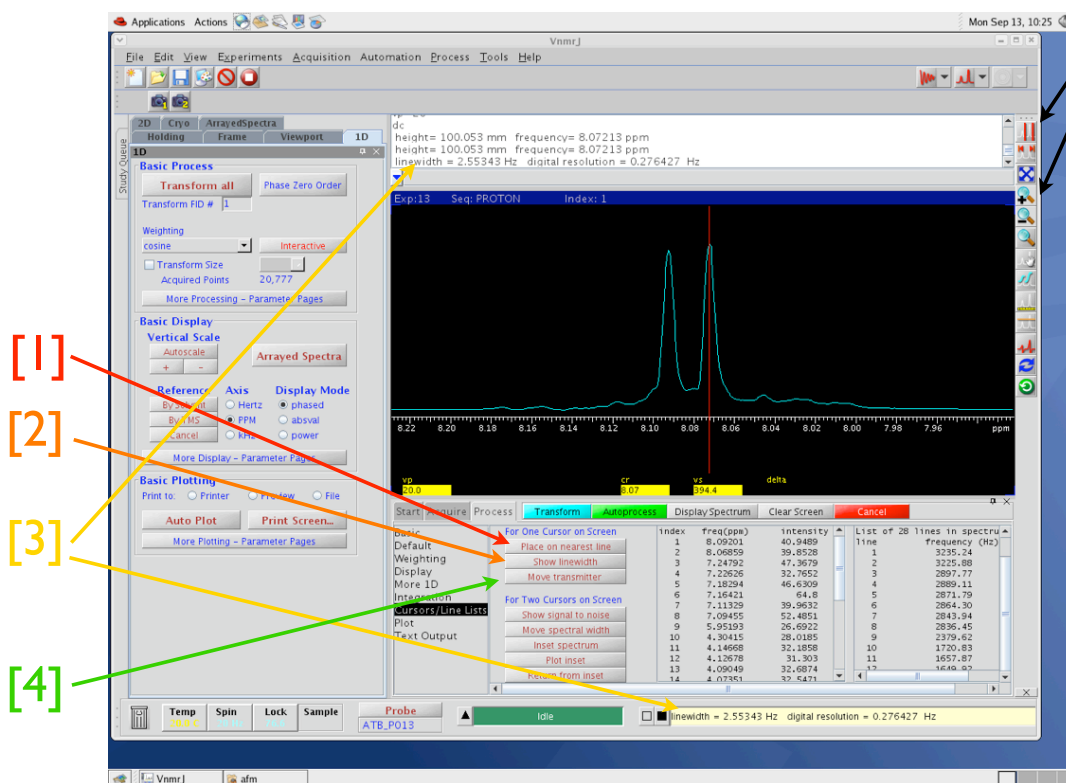


© A.-F. Miller 2010

Figure 8  
Weighting & Linear Prediction

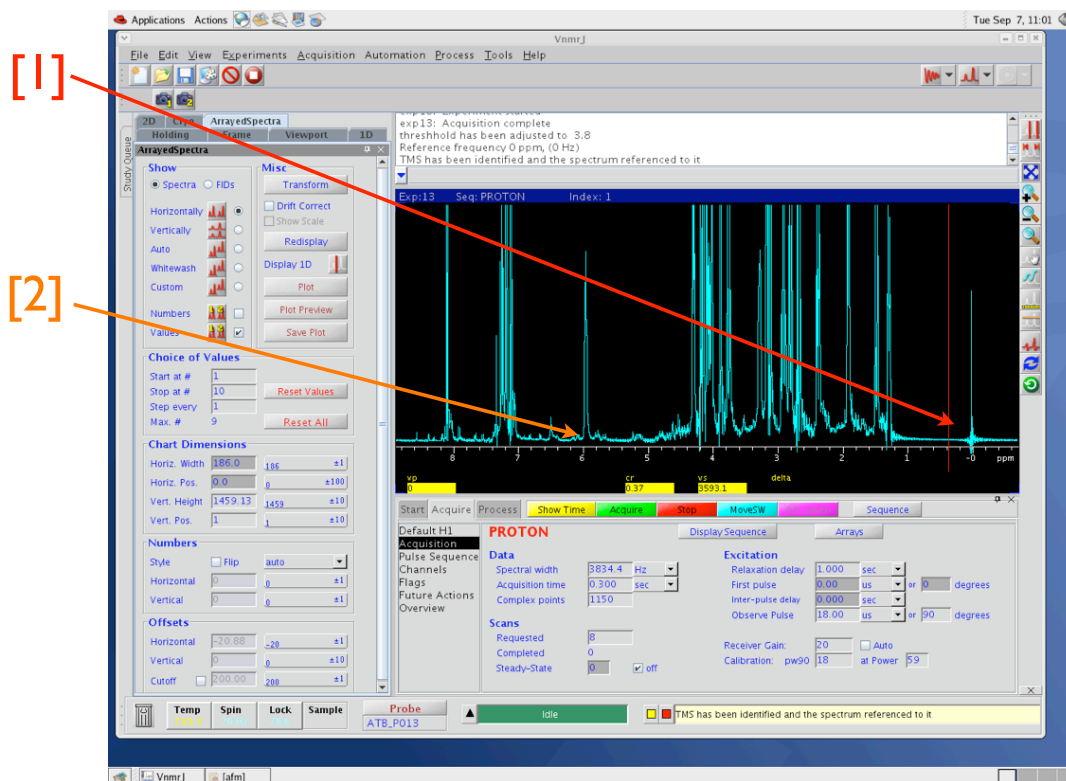


© A.-F. Miller 2010

Figure 9  
Weighting & Linear  
Prediction

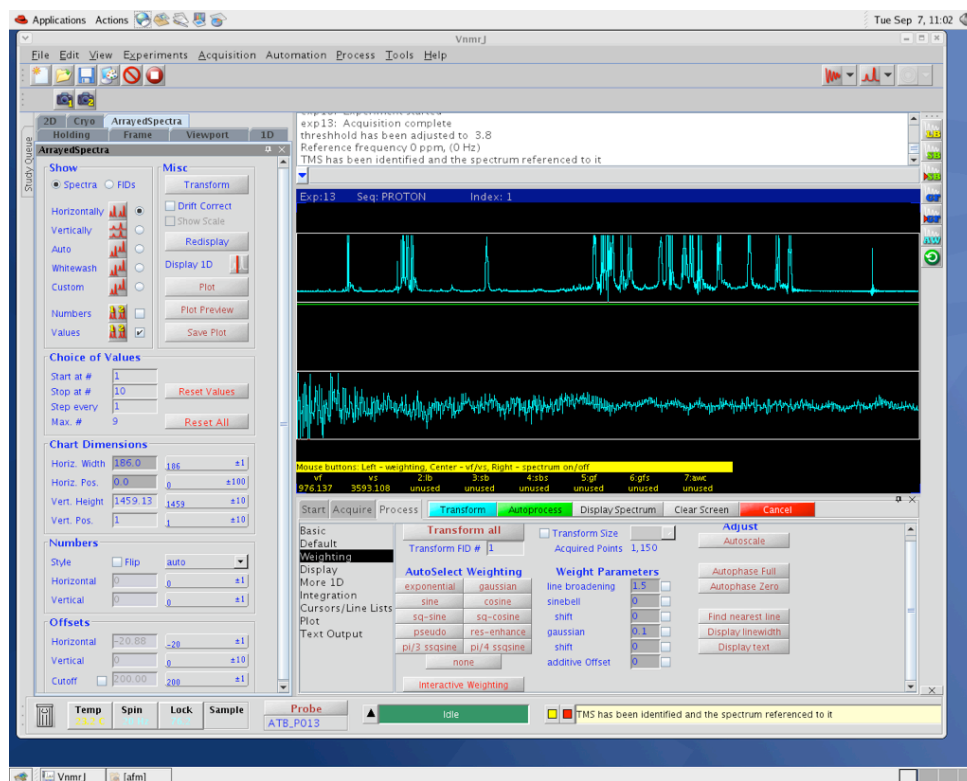

© A.-F. Miller 2010

Figure 10  
Weighting & Linear  
Prediction

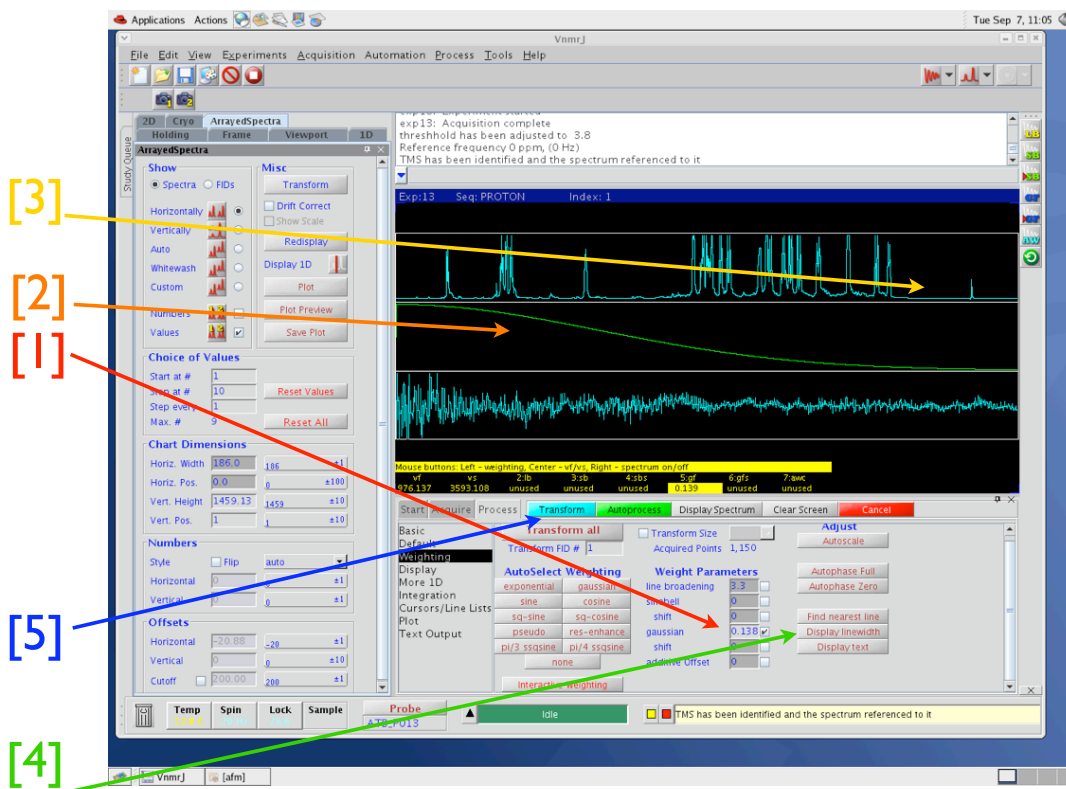


© A.-F. Miller 2010

### Figure 12

#### Weighting & Linear Prediction




Figure 13  
Weighting & Linear  
Prediction

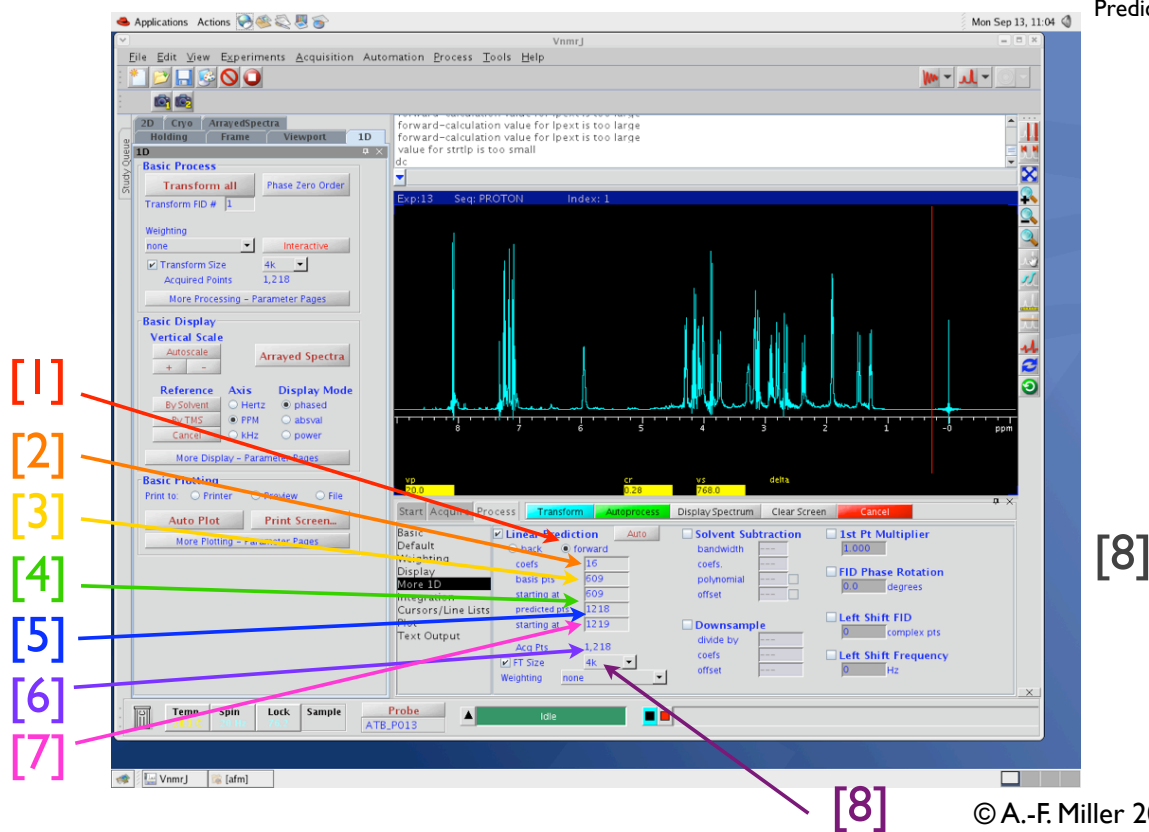


© A.-F. Miller 2010

Figure 14  
Weighting & Linear  
Prediction

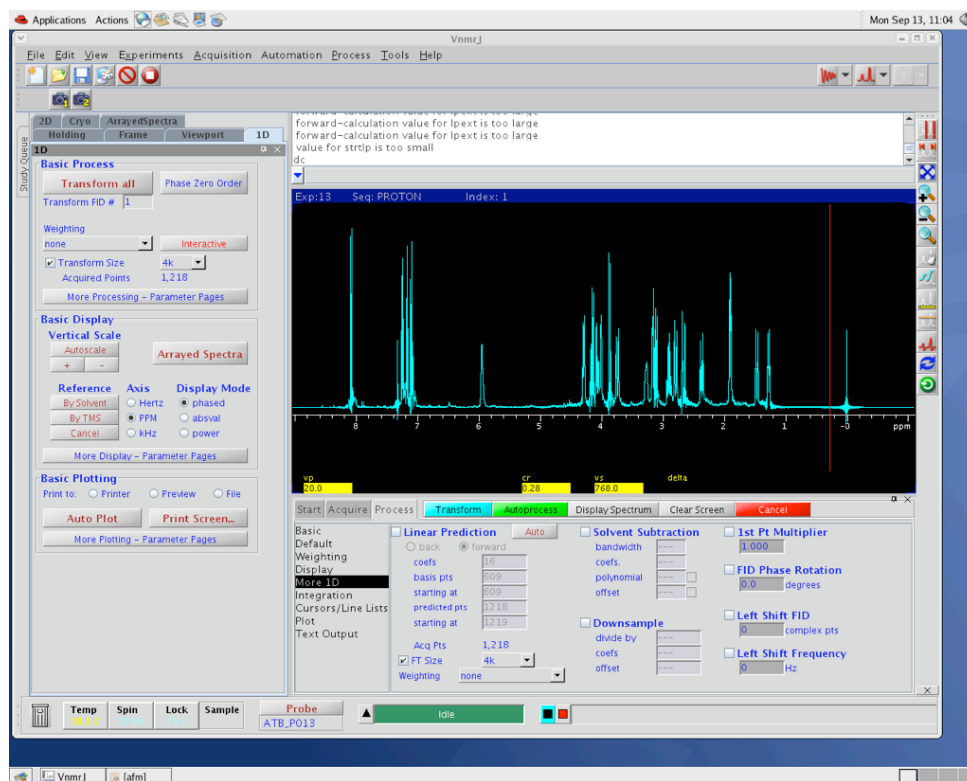


© A.-F. Miller 2010

Figure 15  
Weighting & Linear  
Prediction

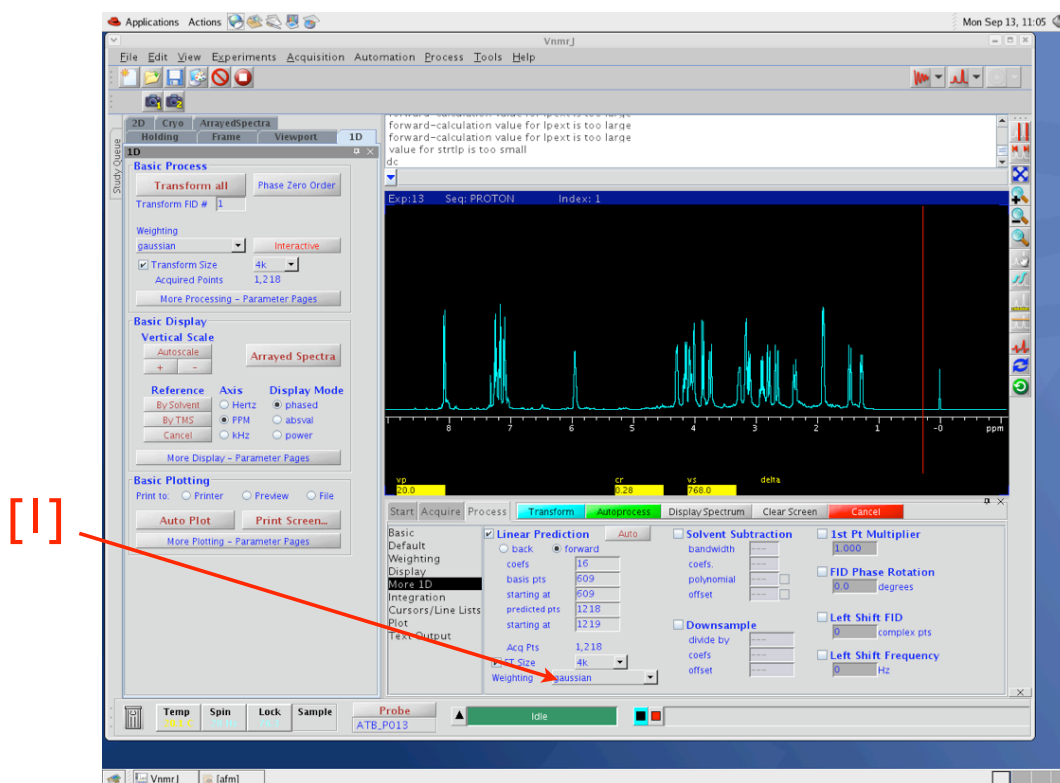

© A.-F. Miller 2010

Figure 16  
Weighting & Linear  
Prediction

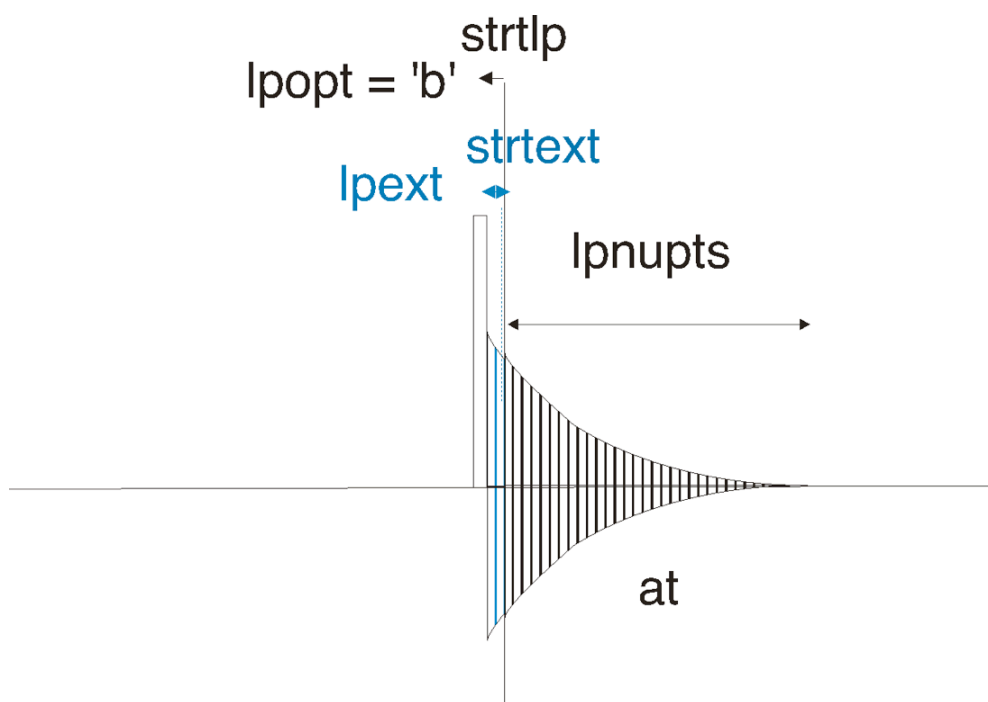

© A.-F. Miller 2010